


SleepStyle™ CPAP SERIES

Ambient
Tracking™
Plus

Specification Sheet

The CPAP range with SmartStick™ Technology

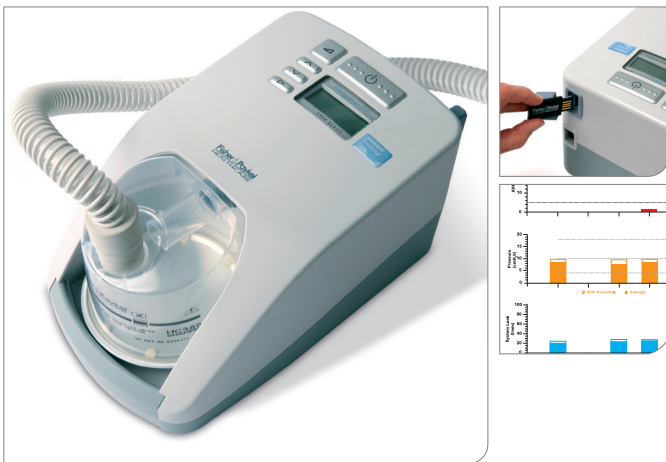

The CPAP range with SmartStick™ Technology offers:

- *Performance monitoring through PerformanceMaximizer™ Software*
- *Two-way communication, including data transfer and remote adjustment of settings, via the SmartStick™*

Fisher & Paykel
HEALTHCARE

The SleepStyle™ 242 and 244, two CPAPs with SmartStick™ Technology.

The **SmartStick™**, utilizing USB technology, provides convenient removable media for access to view and update settings.

This **SmartStick™** Technology, combined with **PerformanceMaximizer™** Software, streamlines therapy analysis with options for compliance data with the **SleepStyle™ 242 CPAP** or full efficacy data with the **SleepStyle™ 244 CPAP**.

PERFORMANCEMAXIMIZER™

Performance monitoring software, **PerformanceMaximizer™**, offers detailed reporting on therapy effectiveness, including compliance, AHI, leak and pressure.


PerformanceMaximizer™:

- is a comprehensive patient management system
- reports any 30 consecutive nights of usage and efficacy using a selectable date range
- reports on compliance of > 4 hours per night, within 30 consecutive days of use
- reports seven nights of detailed efficacy data (SleepStyle™ 244 only)
- reports cumulative data from first use
- allows the printing, saving or emailing of reports
- allows for the easy download of data with the **SmartStick™**

The **SmartStick™** is unique to Fisher & Paykel Healthcare and utilizes USB technology. No additional cables or card readers are necessary.

The **SmartStick™** has 2-way communication, allowing the clinician to view the device settings, compliance and efficacy data, as well as to adjust the device settings.

By utilizing **PerformanceMaximizer™** Software, changes saved to the **SmartStick™** will automatically transfer to the device when the **SmartStick™** is re-inserted into the device by the patient.

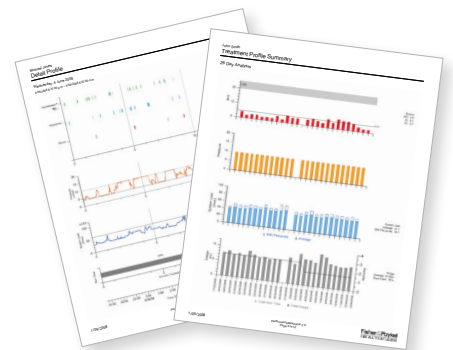



AMBIENT TRACKING™ PLUS

The **SleepStyle™ 200 CPAP Series** is fully integrated with Ambient Tracking™ Plus patented humidification technology. Ambient Tracking™ Plus provides effective humidification under changing conditions by:

- Minimizing condensation to prevent mask pressure fluctuations and algorithm disruptions
- Maximizing humidity during leaks to prevent nasal discomfort



TRUE INTEGRATION

The **SleepStyle™ 200 CPAP Series** is designed for integrated simplicity and enhanced patient comfort. This simplicity ensures quick and easy set-up, use, cleaning and maintenance.

Travel needs are satisfied with a compact and robust CPAP that allows for altitude adjustment and voltage compatibility in most countries of the world.

SleepStyle™ 200 Series CPAP with SmartStick™ Technology				
Performance Category		Performance Feature	242*	244
Pressure Control	A	Altitude Adjustment Control – Manual	●	●
	B	Fully Integrated System	●	●
Humidification Delivery	C	Ambient Tracking™ Plus	●	●
	D	PerformanceMaximizer™ Software	●	●
Performance Monitoring	E	SmartStick™ – Removable Media	●	●
	F	Compliance Data	●	●
	G	Efficacy Data (AHI, Leak, Pressure)		●

**SleepStyle™ 242 model not available in all countries*



SleepStyle™ 200

CPAP SERIES

- **PERFORMANCE** - PerformanceMaximizer™ reports on compliance and efficacy data for efficient patient management.
- **AMBIENT TRACKING™ PLUS** is a unique technology that promotes better humidification and clinical performance superseding conventional humidifiers.
- **CONVENIENCE** - The SmartStick™ allows for access to patient data, utilizing USB technology.
- **EASE OF USE** - Intuitive menu, quiet operation and fully integrated design add up to a quick and easy provider set-up.



The SleepStyle™ 200 CPAP Series will help you provide 24-hour wellness for your patients. Advanced technologies deliver the best in sleep performance for an energized lifestyle. For more information or to view our complete SleepStyle™ CPAP family, visit www.fphcare.com



SLEEPSTYLE™ - 242 AND 244 CPAPS

- HC242** Fully Integrated CPAP with detailed compliance data
- HC244** Fully Integrated CPAP with detailed efficacy data

SLEEPSTYLE™ - 242 AND 244 CPAPS INCLUDE

- CPAP Unit
- Carry Bag
- Water Chamber
- Breathing Tube
- Air Filter
- Water Funnel
- SmartStick™
- SmartStick™ Mailer*

ACCESSORIES

- HC385S** Standard Humidification Chamber
- HC355*** Extended Life Humidification Chamber
- 900HC010** Connector
- 900HC221** Tubing including Connector
- 900HC240** Air Filter (2-Pack)
- 900HC232** Compliance Download Serial Cable
- 900HC611** SmartStick™ (5-Pack)
- 900HC615** PerformanceMaximizer™ Software
- 900HC236** Isolated Serial Adaptor
- 900HC226** Lead Extension Kit
- 900HC225** Water Funnel (10-Pack)
- 900HC630** SmartStick™ Mailer* (25 pack)

PRODUCT SPECIFICATIONS

- Dimensions:** 10.8" x 6.7" x 5.5"
(275 mm x 170 mm x 140 mm)
- Weight:** 4.4 lbs (2 kg)
- Chamber Capacity:** 400 ml
- Electrical Ratings:**
Supply Frequency: 50-60 Hz
Voltage and Current: 100-115 V ~ (1.2 A, 1.3 A max)
220-240 V ~ (0.8 A, 1.8 A max)
- Operating Pressure:** 4-20 cmH2O

* Not available in all countries

FOR MORE INFORMATION, PLEASE CONTACT YOUR LOCAL FISHER & PAYKEL HEALTHCARE REPRESENTATIVE:

Fisher & Paykel Healthcare has a policy of continuous product improvement and reserve the right to alter specifications without notice. Patented in principal countries around the world.

Manufacturer Fisher & Paykel Healthcare Ltd, 15 Maurice Paykel Place, East Tamaki, Auckland 2013, P O Box 14 348, Panmure, Auckland 1741, New Zealand
Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.com Web: www.fphcare.com

Australia Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22 **Benelux** Tel: +31 40 216 3555 Fax: +31 40 216 3554 **China** Tel: +86 20 3205 3486 Fax: +86 20 3205 2132
France Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7181 98599 0 Fax: +49 7181 98599 66 **India** Tel: +91 80 4284 4000 Fax: +91 80 4123 6044
Irish Republic Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205 Fax: +81 3 3661 7206 **Northern Ireland** Tel: 0800 132 189
Spain Tel: +34 902 013 346 Fax: +34 902 013 379 **Sweden** Tel: +46 8 564 76 680 Fax: +46 8 36 63 10 **Switzerland** Tel: 0800 83 47 63 Fax: 0800 83 47 54
Taiwan Tel: +886 2 8751 1739 Fax: +886 2 8751 5625 **UK** (EU Authorized Representative) Fisher & Paykel Healthcare Ltd, Unit 16, Cordwallis Park, Clivemont Road, Maidenhead, Berkshire SL6 7BU, UK Tel: +44 1628 626 136 Fax: +44 1628 626 146 **USA** Tel: 1 800 446 3908 or +1 949 453 4000 Fax: +1 949 453 4001